### **Technical Specification**

**Mode Type** 

Whole Blood, Predilute, Capillary

**Total Parameters: 22** 

CBC: WBC, RBC, HGB, HCT, MCV, MCH, MCHC, RDW-CV, RDW-SD, PLT, MPV, PCT, PDW-CV, PDW-SD, P-LCR, P-LCC.

3 Part Diff: Lymph % & #, Mid % & #,

Gran% & #

Principle of measurement

RBC/PLT/WBC/ DIFF: Electrical Impedance Hb: Cyanide Free Colorimetry

*MCV*: Measured *HCT*: Calculated

Histograms
3 (WBC/RBC/PLT)

Sample Volume

Whole Blood : 9 µL Pre-diluted : 20 µl Capillary : 20 µl **Linearity Range** 

WBC :  $0-300 \text{ (x } 10^3/\mu\text{L)}$ RBC :  $0.00-8.50 \text{ (x } 10^6/\mu\text{L)}$ Hb : 0-25.0 (g/dL)HCT : 0-67 (%)PLT :  $0-3000 \text{ (x } 10^3/\mu\text{L)}$ 

Calibrator

ELite H Cal (3mL)

Open Vial Stability at 2-8 °C: 7 Days

**Tri-level Controls** 

Elite H3 CON L, N, H (3mL) Open Vial Stability at 2-8 °C: 16 Days

Thoroughput

60 Tests/Hr

**Quality Control** 

500 files with L-J graph X-Bar Analysis

**Data Storage** 

50,000 Results with Graphs

Dimension (mm) 364x477x417

Weight (Kg)

25

Reagents

Erba H 360 Dil ERBA H 360 Lyse Elite H Clean

Interfaces

4 USB + 1 LAN Port

**Operating Environment** 

Temperature : 15-35 °C Humidity : 30-85%

**Power Requirement** 

A.C.100-240V 50/60Hz 200VA

### Transasia - Erba Group

# Global footprint



### **India Advantage**



Installations



Customer Support
Least downtime





Supply Chain Widest Network in the Industry

Sales & Marketing Largest in the IVD Industry

Manufactured by:

Erba Lachema s.r.o. Karásek 2219/1d, 621 00 Brno, Czech Republic Tel.: +420 517 077 111

E-mail: diagnostics@erbamannheim.com

1800 130 8226
W: www.erbamannheim.com
Stay connected: in 

■ C

Toll Free Number

### TRANSASIA BIO-MEDICALS LTD.

Transasia House, 8 Chandivali Studio Road, Andheri (E), Mumbai - 400 072. T: (022) 4030 9000, F: (022) 2857 3030 E: responses@transasia.co.in









A product of Erba Europe

### **CONVENIENT**


360

- 10.4" large coloured LCD touch screen with graphic icons
- Single page view of full report
- Convenient daily maintenance with single shut down procedure



## PRECISE & ACCURATE

- Excellent short & long term precision
- Highly accurate with results correlating with microscopic examination & clinical findings
- Specific flagging messages to alarm for microscopic slide review

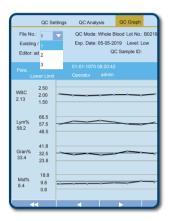

### **MORE WITH LESS**

- Just 9 µL sample volume: all age groups
- P-LCR, P-LCC, PDW-SD, PDW-CV aids differential diagnosis of cases with abnormal platelet counts
- RDW-SD, RDW-CV & MCV aids in screening for causes of anaemia



# COMPACT AUTOMATED HEMATOLOGY ANALYSER WITH 22 REPORTABLE PARAMETERS



RIMS

Real time Inventory

Management System with RFID





Erba

Transcribe other test results - consolidate reports

### **ADVANCED SOFTWARE**

- Large storage capacity: 50,000 reports with graphs
- 4 USB ports Bar code scanner, external printer, key board, pen drive
- 1 LAN port for bi-directional interfacing (Optional LIS solution)

